




# ESK-001/6

Hydraulic and manual bed

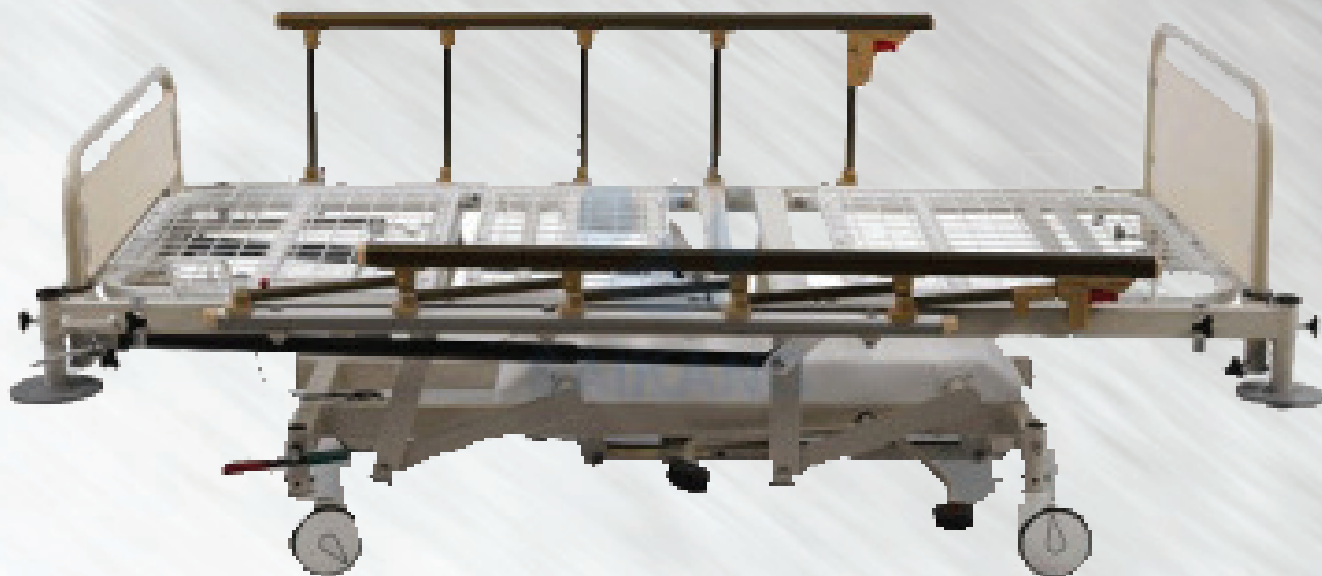

MEDICAL FURNITURE



# ESK-001/6

Hydraulic and manual bed

## Description:

This is a five-function manual and hydraulic bed. Turn the right crank, the backrest will be up and down, turn the left crank, the knee part will be up and down. By pressing the UP/DOWN pedals, bed height is adjustable. Trendelenburg & Reverse-trendelenburg function are also available.

- Detachable head/foot board, round bumpers at four corners of the bed.
- Aluminum alloy side rail, 5-staff, collapsed type with safe lock.
- Bed sheet top made from mesh steel, epoxy powder coated and baked finish.
- Bed framework made from mild steel ERW rectangular tubes.
- Central locking foot pedal control with casters diameter 125mm.
- With mattress

## Technical Configuration:

Overall size: 2160\*990\*440-845mm  
 Height adjustment (without mattress): 440-845mm  
 Backrest lifting angle: 0-80°  
 Kneerest lifting angle: 0-45°  
 Trendelenburg: 25°  
 Reversed trendelenburg: 25°



ESSE3 srl, Via Garibaldi 30  
 14022 Castelnuovo D.B. (AT)  
 Tel +39 011 99 27 706  
 Fax +39 011 99 27 506  
 e-mail [esse3@chierinet.it](mailto:esse3@chierinet.it)  
 web: [www.esse3.dreamgest.com](http://www.esse3.dreamgest.com)

